

Image instructions



funke-medical.de/image-instruction/cubeseat-en.pdf




CUBE SEAT

Anti-decubitus cube cushion

- Autoclavable according to RKI*¹ guidelines
- Applicable up to stage 2 (EPUAP)
- Middle risk (Norton-,Braden-,Waterlow-Scale)



The CUBESEAT is an anti-decubitus cube cushion that specifically helps to relieve pressure. Thanks to its innovative cube structure, seating comfort can be individually adjusted. The foam cubes can be replaced with visco foam cubes as required to customise seating comfort.

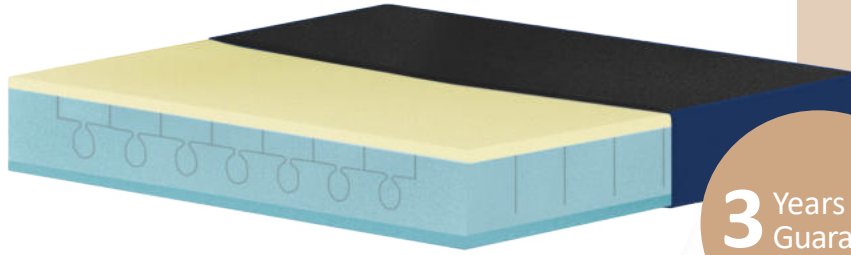

Characteristics mesh cover

- Removeable, washable up to 60°C
- Mesh³ fabric guarantees an excellent microclimate

Available with PU cover

- Removeable, washable up to 95°C, disinfection by wiping, disinfectant-proof*³

Special Features

- Seating comfort individually adjustable
- Unique cube structure for individual adjustment of seating comfort
- Innovative contour ensures excellent flexibility and stability of the foam cubes
- Non-skid bottom

*1 Robert-Koch-Institute (RKI) is a part of the Federal Ministry of Health in Germany. It is responsible for disease prevention and control.


*2 Guarantee applies to functional properties of the cushion core (pro rata basis).

*3 For more details see instructions for use.

Size

(Width x Length x Height)

Size	Item No.
38 x 40 x 8 cm	C-SEATM3840
40 x 38 x 8 cm	C-SEATM4038
40 x 40 x 8 cm	C-SEATM4040
40 x 42 x 8 cm	C-SEATM4042
42 x 40 x 8 cm	C-SEATM4240
42 x 45 x 8 cm	C-SEATM4245
45 x 42 x 8 cm	C-SEATM4542
44 x 45 x 8 cm	C-SEATM4445
45 x 44 x 8 cm	C-SEATM4544
45 x 46 x 8 cm	C-SEATM4546
46 x 45 x 8 cm	C-SEATM4645
42 x 50 x 8 cm	C-SEATM4250
50 x 42 x 8 cm	C-SEATM5042

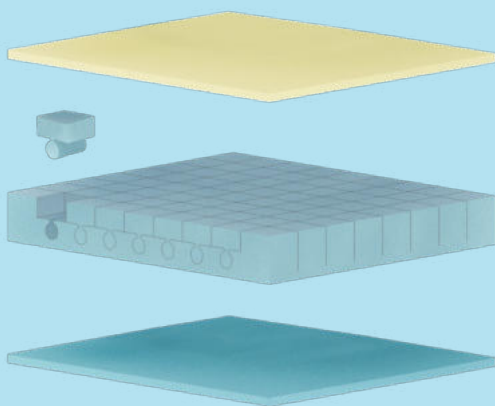

Structure

- Top layer made of visco-foam for an ideal adaption to the contours of the body
- Middle layer consists of unique cube structure for individual adjustment of seating comfort
- HR foam bottom layer for increased stability
- Inner pocket for top plate/bottom plate
- Visco-foam cubes available as accessories

Technical data

- Breathable open-cell foam
- Autoclavable
- Waterproof but breathable
- Individually adjustable
- Medical device according to MDR
- Maintenance-free
- Patient weight: 0-120 kg

Size foam cube
(Width x Length x Height)
5 x 5 x 4,2 cm

